




Installation description –

How to mount a 2- or 4-point belt to a chair having a dynamic seat plate

Item Kit = 45278 consists of		
Amount	Item number	Description
2	82794	Hip belt bracket for seat plate Dynamic
2	20296	Square nut M6
2	29089	Hexagon bolt socket flat head M6 X 12
4	20063	Washer M6
4	20058	Lock nut M6
4	20800	Bolt Hex socket button Head M6 X 16

At a chair with a dynamic seat plate, you will find already item 82794 installed.
If you miss the hip belt bracket, please order the kit 45278 twice.

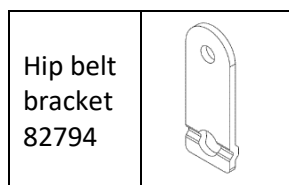

To add the hip belt bracket to the dynamic seat plate, add the square nut 20296 to the rail of the seat plate and fix the bracket from outside with the Hexagon bolt 29089.

Together with the 2- or 4-Point belt you will get the below shown Flat-Mount brackets.
Use the remaining parts to mount the Flat-Mount bracket to the Hip belt bracket.



Please follow the instruction on the following side to adjust the belt.

Two-Point Pelvic Support Belt

1. Locate desired anchor point.
2. Center the belt on the pelvis and thread the webbing end through the Flat-Mount bracket.
3. Cinch the webbing snug.



Four-Point Pelvic Support Belt

1. Locate desired anchor point
2. Follow the installation instructions above to fit the primary straps.
3. Attach the secondary straps to the chair so that they anchor the primary straps.
4. Tighten secondary straps to prevent movement of the primary straps.



ADJUSTMENT:

When properly adjusted and the belt tightened, it should fit snug so that the user's pelvis is secure. On padded belts, tighten the D-ring pull until all slack is removed and the pads under the buckle is touching. If length adjustment is necessary to achieve a good fit or to remove the belt for cleaning, this can be done using the belt end fittings.



SAFETY CHECK:

When properly positioned in a wheelchair, have user lean forward and side to side to check fit.

Check for:

1. Normal operation of buckle and adjustment straps.
2. Comfort - Look for areas of irritation.
3. Position -If too high or too low. Adjust anchor points.
4. Interference with other devices - Relocate anchor points as necessary.

This instruction is part of every 2- or 4-point belt and must, therefore, be forwarded to the user or carer and should be stored.

If you have any further question related to this item or mounting description, please contact our customer service desk.

These and other mounting description can also be found on our homepage at My-netti.com.